


Calibra 2 PRO P



Suggested fitting range:

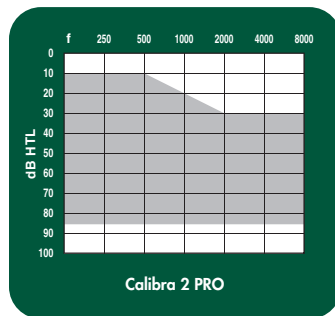

Calibra 2 PRO P BTE

Technical Description:

The Calibra 2 PRO P (power) directional BTE is a premier-level 100% DSP circuit with seven channels, four adjustable compression channels, three independently programmable memories with a fourth memory dedicated solely for audio input. This BTE features an endless volume control wheel. Other outstanding features of the Calibra 2 PRO P include Automatic Microphone Management, Broadband Adaptive Directional Microphone, DataLogging, Feedback Cancellation, and four personalized choices for Environment and Noise Management: speech, speech in noise, noise, and music. Rexton's easy-to-use software helps you program this superior technology by incorporating new and improved fitting protocols in a familiar programming platform.

Fitting Recommendation:

The Calibra 2 PRO P BTE is an excellent choice for mild to severe hearing loss of many configurations including flat, sloping, high frequency, and reverse curve.

Standard Features:

- 7 Channels
- 4 Memories, one dedicated for audio input
- Volume Control
- Automatic Adaptive Directional Microphone
- Feedback Cancellation
- Personalized Noise Management
- DagaLogging
- 4 AGC Adjustments
- Power On Delay
- Memory Selection Tones
- Push Button On/Off
- Telecoil
- Compatible with wireless products, cell-phones, FM systems, and MLX
- Beige, grey, brown, or metallic champagne housing colors
- 13 A battery type housing

Available Options:

- Eyeglass adapter
- Small Earhook
- Audio Shoe/Integrated Battery Door

ANSI S3.22-2003

MAXIMUM POWER OUTPUT	Peak	dB OSPL90	130
	HF-average	dB OSPL90	126
GAIN (Input at 50 dB)	Peak	dB	65
	HF-average	dB	58
REFERENCE TEST GAIN		dB	48
FREQUENCY RESPONSE		Hz	100-6600
COMPRESSION	Attack time	ms	VARIABLE
	Release time	ms	VARIABLE
TOTAL HARMONIC DISTORTION	Maximum	%	7
EQUIVALENT NOISE LEVEL		dB SPL	16
INDUCTION COIL SENSITIVITY	HFA-SPLITS@31.6 mA/m		
	Left/Right	dB SPL	105/109
BATTERY DRAIN	at RTG., Input 65 dB SPL	mA	0.9
BATTERY LIFE	Type 13 zinc air	hrs	~250

Rexton, Inc.

5010 Cheshire Ln. N., Suite 2
Plymouth, MN 55446

800-876-1141

Fax: 763-553-9129

www.rexton-online.com

REXTON

Calibra 2 PRO Power BTE

